



# EASYCAP

EEG Recording Caps and Related Products

EASYCAP GmbH  
Steingrabenstrasse 14  
DE-82211 Herrsching  
Germany  
Delivery Address: Am Anger 5, DE-82237 Woerthsee-Etterschlag

Tel +49 (0) 8153 88702-00  
Fax +49 (0) 8153 88702-10  
www.easycap.de  
info@easycap.de

## 96Ch Standard Cap for actiCAP for actiCHamp


Green Holders: Label 1-32, Hardware Channel 1 – 32

Yellow Holders: Label 1-32, Hardware Channel 33 - 64

Red Holders: Label 1-32, Hardware Channel 65 - 96

Black holder: Label & hard-wired Gnd

### Electrode Names and Number Labels

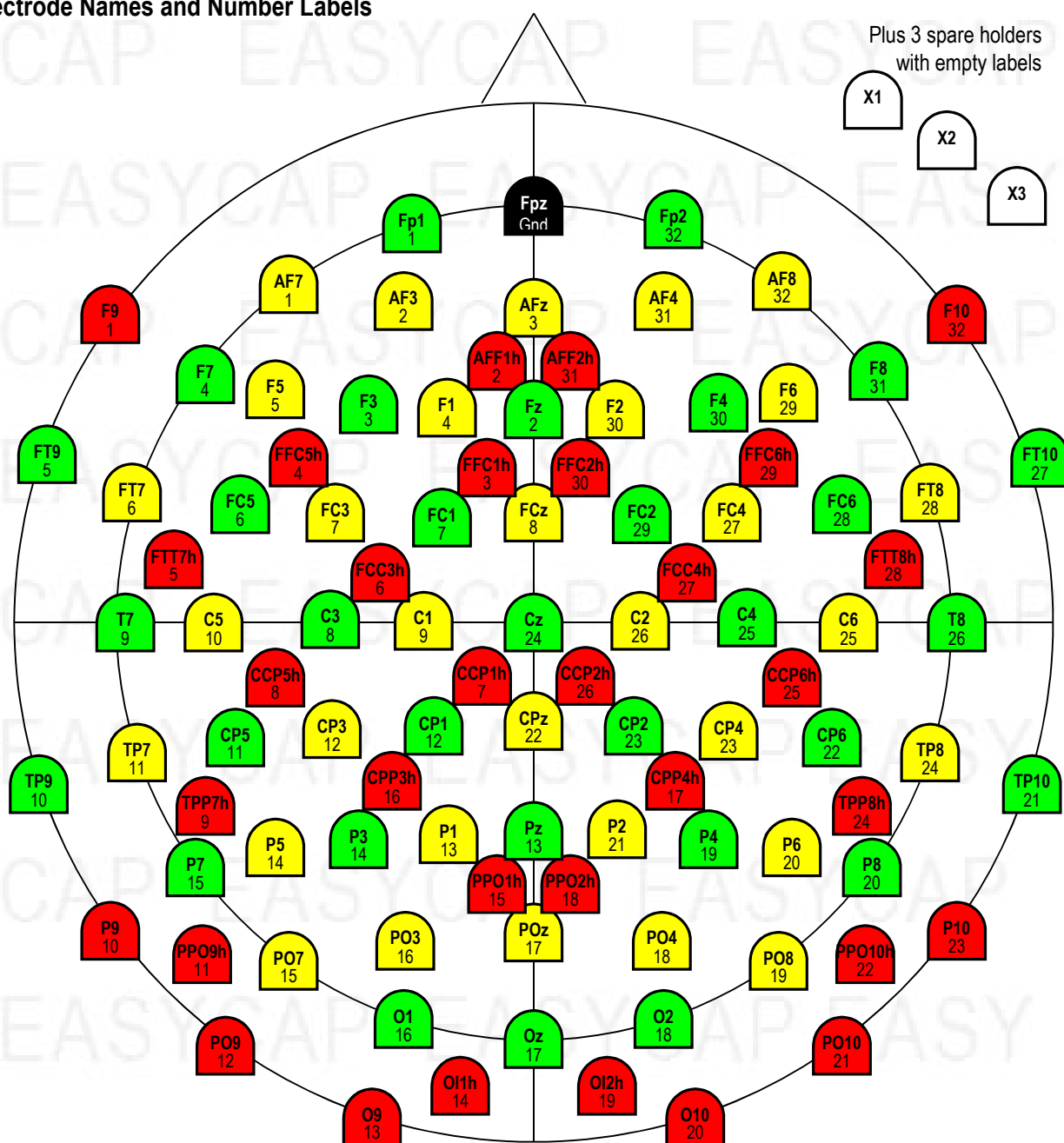

Y:\AL\Hauben\_Layouts\AC\_AP\_actiCAP\_actiCHamp\APAP-96VAP-96.docx

See [actiCHamp 160Ch guided Standard.doc](#) for a guided map of the numerical channel number order.



EasyCap GmbH, Steingrabenstrasse 14, DE-82211 Herrsching, Germany

## Details

### Ordering Information

For ordering please give **Article Number, Cap Cut, and Size**

(Example: CAP-96, Caucasian, 58):

- Article Number: **CAP-96**
- Cap Cut: **Caucasian** or **Asian**
- Size (given in cm head circumference):
  - Adult caps: **54, 56, 58, 60, 62, 64** (average male: 58, average female: 56)
  - Children caps: **52** (5-10 years), **54** (11-14 years)

The catalogue-number comprises the cap as described, serial number, and this document; all packed in a labelled cardboard box. For further information about accessories or consumables, please visit our website or contact our local distributor.

### Cap

White Subtemporal Cap with integrated chin belt

Sizes 52-64 made from High Precision Fabric

Options: *Caucasian* or *Asian*, *Size*

The labelled holders are mounted according to the figure (1-96 transparent, Gnd black). The positions are chosen from 10%-System and from 5%-System<sup>1</sup>. Positions TP9/TP10 approximately meet left and right mastoids.

Each cap comes with 3 additional holders (apricot empty label) for use with double-sided adhesive rings to place electrodes on bare skin. The grip of the washer is better, when both, skin and holder, are degreased with alcohol prior to application.

### Accessories

The cap can be fixed either by the integrated chin belt alone, or additionally or alternatively by a chest belt (not included). Normally the decision between these 2 methods can be made according to what is more comfortable for the subject. Under difficult conditions the chest belt will hold the cap in place a little better.

### Service

In case of any further questions please contact our local distributor or [info@easycap.de](mailto:info@easycap.de). Please have the serial number of your cap at hand (syntax Axxxxx, located on the backside label of the cap).

---

<sup>1</sup> Electrode Nomenclature according to:

Oostenveld, R. & Praamstra, P. The five percent electrode system for high-resolution EEG and ERP measurements. *Clinical Neurophysiology* 2001; 112: 713-719