



### Support Tip

## Have you located the Solutions help documentation for Analyzer 2?

by Dr. Roland Csuhaj

The Solutions are small programs which add useful functionalities to BrainVision Analyzer. For example, one of them enables you to import files saved in the BESA file format, while another makes it possible to export the area under the curve around the peak markers.

Since the Solutions are not part of the main BrainVision Analyzer 2 program, they are not described in the User Manual. Each Solutions file contains its own associated help documentation. In Analyzer 2, the Solutions / Help / Solutions Help menu opens the Solutions Help Explorer program,

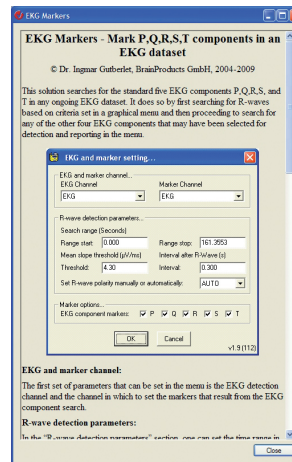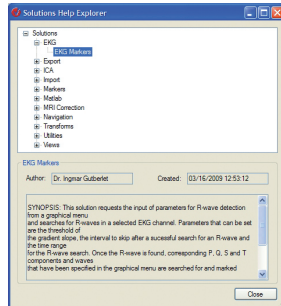

where a short description appears when you left-click on an item.

The short descriptions generally consist of just a few sentences, while double-clicking opens a separate window that displays more detailed help documentation.

Here is where you can download the latest version of the Solutions for BrainVision Analyzer 2:  
[www.brainproducts.com/downloads.php?kid=9&tab=3](http://www.brainproducts.com/downloads.php?kid=9&tab=3)